

CHARGING INSTRUCTIONS

1. Plug USB charging cord into port connector on the device.
2. Plug USB charging cord into a USB power adapter, and into outlet for charging.
3. Charge your Illumina™ for 2-4 hours.
4. Device is fully charged when indicator light is solid green.

LIGHT SENSITIVITY TEST

Test yourself for light sensitivity prior to initial treatment:

1. Power on your device.
2. Place device on or within 1/4" of the skin on your forearm for 3 minutes.
3. If area has turned red and lasts for more than 2 hours your skin is light sensitive. Discontinue use and contact **Customer Service** at **877-425-4468** or email **support@staraligners.com**.

For more information, please review the **Illumina™ User Guide** provided with your device.

Manufactured exclusively for StarAligners

9920 Pacific Heights Blvd, Suite 150
San Diego, CA 92121

ILLUMINA™
ORAL THERAPY SYSTEM
BY STARALIGNERS

Made in China

(877)425-4468 | staraligners.com | [#straighttohappiness](https://twitter.com/straighttohappiness)

PATIENT GUIDE

ILLUMINA™
ORAL THERAPY SYSTEM

RELIEVES TEMPORARY ORAL DISCOMFORT

TEETH WHITENING ACCELERANT

ENHANCES GUM HEALTH



BY STARALIGNERS

ILLM022PB-R16001

Getting Started

Get the most out of your oral care program!

The StarAligners Illumina™ Oral Therapy System uses red and infrared light to provide temporary relief from oral discomfort. Using low-level blue light, it whitens teeth when used in conjunction with peroxide gels or foams.

SYSTEM INCLUDES

- Illumina™ Oral Therapy Device
- Power Adapter
- USB Charging Cord
- Storage Bag
- User Guide
- StarAligners Teeth Whitener + Aligner Cleaner



HOW IT WORKS

The StarAligners Illumina™ Oral Therapy System uses red and infrared light to provide temporary relief from oral discomfort. Using low-level blue light, it whitens teeth when used in conjunction with peroxide gels or foams.

Relieves temporary oral discomfort

Teeth whitening accelerant

Slows and in some cases, reverses receding gumlines

Provides increased oxygen and nutrients to the gums

Reduces the effects of gingivitis


Promotes healing of gum tissues

OPERATING INSTRUCTIONS

1. Fully charge the Illumina™ device for 2-4 hours.
2. Press and hold the On/Off button for 2 seconds to see lights appear.
3. Place silicone tray into mouth, securing device into place.
4. Use device for 10 minutes for gum health and teeth whitening.
Note: The device will auto-shut-off after completion of 10-minute self-treatment.
5. Remove and repeat the opposite arch.
6. Device may be used multiple times per day, up to 5 treatments per day.
7. Place Illumina™ device into protective storage case when not in use.

TEETH WHITENING

- Self-treatment time is 10 minutes per arch per day.
- Device will automatically shut off after 10 minutes.
- Use with StarAligners teeth whitening foam to enhance teeth whitening.

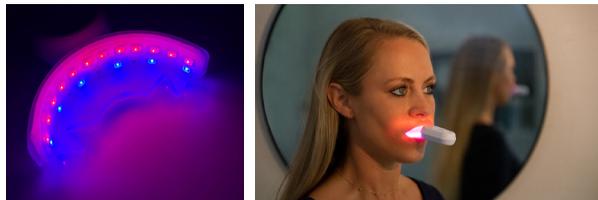

GUM HEALTH

Use once daily for minimum of 10 minutes per arch.

MAINTENANCE

1. This is a medical device. It is important to keep it clean.
2. Charge device only with supplied USB power adapter.
3. Power indicator light will begin to flash when battery levels are low, indicating it is appropriate to begin charging.

CLEANING

1. Turn you device off and unplug from any power source.
2. Gently wash mouthpiece in lukewarm water or Aligner Cleaning Foam.
3. Let device air dry.